

ScopeClip

New single-use Hemoclip

Different advances have been made to this latest generation of haemostatic clips, with the following objectives:


- ▶ A better grip on the tissues.
- ▶ A greater closing force.
- ▶ The possibility of extraction, if necessary




Repositionable



Rotatable

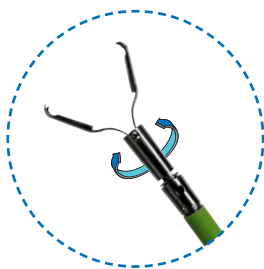

Removable

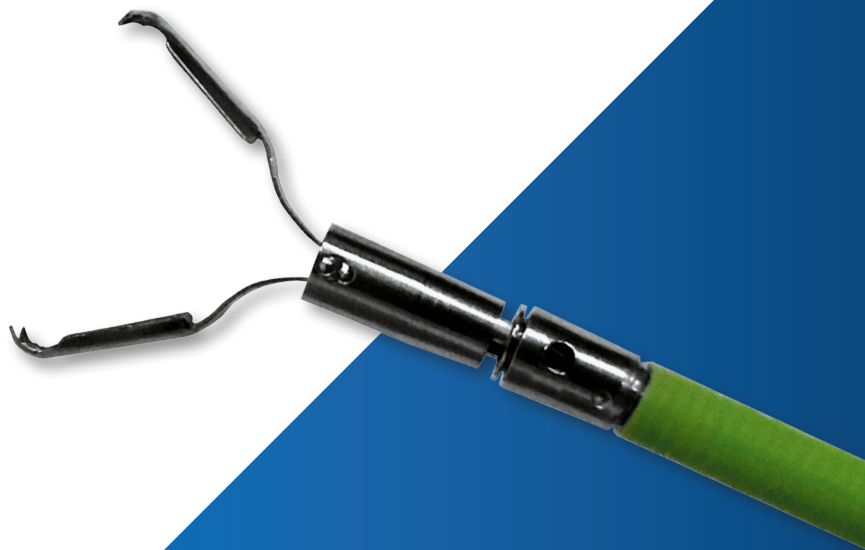


ScopeClip

New single-use Hemoclip


ScopeClip - Hemoclip rotative, repositionable and removable single-use

Reference	Clip opening (mm)	Clip angle	Catheter Ø (mm)	Useful length (mm)	Color code	PKG
HC-23091355108	9	135°	2,6	2300	Dark blue	5
HC-23111355104	11	135°	2,6	2300	Green	5
HC-231313551018	13	135°	2,6	2300	Violet	5
HC-231613551044	16	135°	2,6	2300	Ligth blue	5

Characteristics

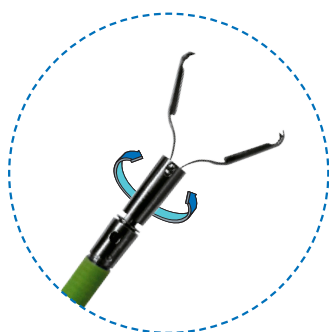

▶ **Rotatable:** 360° rotation in two directions for a better angle of subjection offering a greater manageability and an exact positioning.

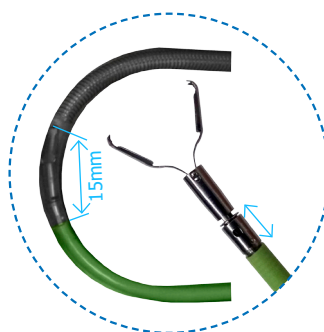

▶ **Shorter distal end:** to facilitate the progression of the clip through the working channel when the endoscope is strongly supported.



▶ **Repositionable:** Possibility to open and close the clamp up to 5 times for a more precise positioning.



▶ **Better tissue grip and greater closing force:** Prevent tissue leakage when closing the clamp and increase the degree of tissue compression for better vascular closure.



▶ **Removable:** If necessary (the clip is incorrectly positioned, the clip does not easily detach from the tissues at a distance from its placement), the clip can be reopened and removed with a polypectomy handle, without damaging the tissues.

Description: Single-use haemostatic clip. **Indications:** Endoscopic placement of clips for haemostasis or marking. **Class:** IIb. **Notified body:** CE0123. **Manufacturer:** Hangzhou AGS MedTech Co, Ltd. **Proper use:** Before use, carefully read the instructions in the package insert or on the device label. **Version:** SCL_E_5. Rev.05. Date 15-07-2022

Creo Medical authorised local distributor

Creo Medical UK Ltd

Chestnut House, N.Minster Business Pk.

YO26 6QR York, United Kingdom

Email: info.uk@creomedical.com

Tel: +44 (0) 1904 786888

Technical changes and mistakes reserved. Images may show non standard options.

